


CE 0197

MEDIKABAZAAR

MANI®

FOR ADVANCED TREATMENT

PEESO REAMERS GATES DRILLS

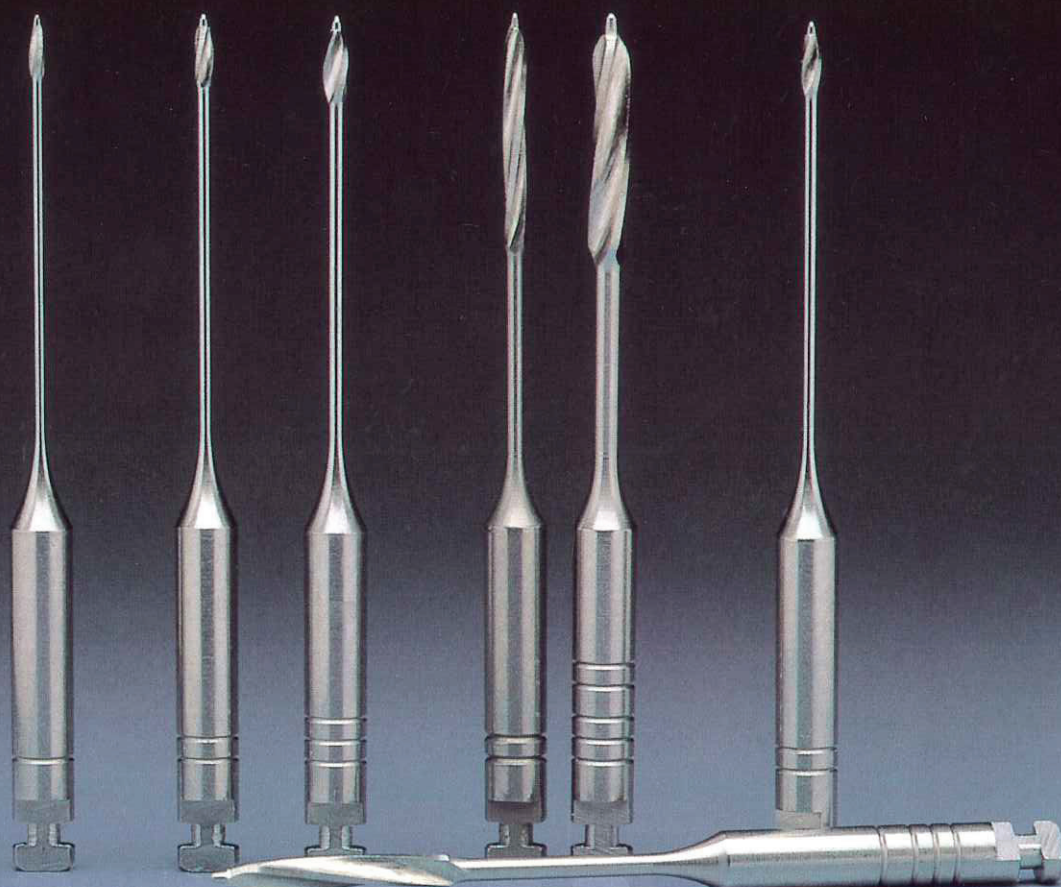

MANI, INC.



For advanced treatment ——— We serve with micro processing technology

PEESO REAMERS GATES DRILLS

**Get
rusted?**

The material developed by MANI has made it possible to commercialize safe and sharp **PEESO REAMERS & GATES DRILLS** with corrosion resistance to any kind of sterilization including autoclave, chemical vapor, and super acid water.

The anti-rust material

We adopted "MANI hard-fiber stainless steel" which has uniform and stabilized hardness as the raw material of this instrument.

MANI Hard-Fiber Stainless Steel
(Austenitic Stainless Steel)



Martensitic Stainless Steel

Length : 28mm, 32mm and 38mm

In addition to standard 32mm length, 28mm and 38mm lengths are available.

You can choose the length exactly you need for molar or canine etc.

PEESO REAMERS

size	φmm
1	0.70
2	0.90
3	1.10
4	1.30
5	1.50
6	1.70

Length(L) : 28,32,38mm
Packed : 6pcs. to a package
Size : 1,2,3,4,5,6 Assorted 1~6

GATES DRILLS

size	φmm
1	0.50
2	0.70
3	0.90
4	1.10
5	1.30
6	1.50

Length(L) : 28,32,38mm
Packed : 6pcs. to a package
Size : 1,2,3,4,5,6 Assorted 1~6

Safety design

Safety guide of PEESO REAMERS and GATES DRILLS prevent ledges or perforations. Designed to be broken at very near to shank.

Easy identification

The size of PEESO REAMERS and GATES DRILLS can be easily identified with the number of grooves on the shank.



CAUTION

Before using the MANI PEESO REAMERS and MANI GATES DRILLS :

●For dental use only. ●Sterilize PEESO REAMERS and GATES DRILLS before use. ●Insert the shank into the handpiece according to the instruction of the handpiece manufacturer. The insertion should be to full shank length. ●This product should be treated as medical waste when disposed.