

CHAMUNDI Left to Right
DENTAL CHAIR MOUNT UNIT




Chamundi Dental Unit Left & Right Hand Operating :

- **Programmable Dental Chair** with 6 programme, 4 programme erasable, gargle, zero.
- **5" TFT based display** having 480/800 resolution, control panel screen with option for 6 programmable chair position.
- **Memory Foam:** Seat, backrest and head rest are fitted with Most latest memory foam covered with Ultra Leather from USA.
- **Light:** Unit mount, LED ("FARO" MAIA) Operating light having On/Off & Intensity control by (un-touch) sensor.
- **Head Rest:** Having lateral left to right movement, rising and lowering, rotation on central axis with pneumatic rigid lock.
- **Hand Rest:** Euro moulded rigid, Dropping down for easy access of patient.
- **Spittoon:** Full Attachment moveable from left to right without using any tool & totally free movement, complete stable.
- **Spittoon Component:** Spittoon move automatically over the body, synchronised with gargle, and goes back to home position when the chair goes to working position. Also can be moved manually.
- Oval shaped bowl to cover more area for easy spitting, detachable spittoon, autoclavable, Special shape water component, Most graceful and clean.
- **Modular System:** Modular arm front loading moveable by 270 degree left to right with mounted on bottom of the seat.
- **Trolley Arm:** Heavy duty arm hinge with manual lock, well balanced.
- Over hanging delivery system, perfect locking by pull, stable, unlocking by pulling and leave the tubing.
- All Valves for air and water are solenoid controlled, Electrically operated 24 VDC valve for life long service.
- **NSK Brushless Optic Micromotor** - NLX Nano, fiber optic with internal spray having chipair function, supplied with NSK 1:5 speed increase contrangle handpiece having optic light and internal spray, push button type.
- **Airotor Points - 2Nos :**
 - a). Super torque Optic **NSK- Dyna LED** handpiece - 1 No
 - b). Mini head handpiece **NSK-Pico** smallest head in the world - 1 No
- 6 way syringe warm water and air - 1No. (LUZZANI Italy)
- 3 way syringe - 1No
- NSK Optic Scaler with 3 scaling tips - 1 Set
- **Multi-Function Foot Control.**
- **Wetline Cuction:** continuous rating 300 L/Minute - 1No
- **Assistant Console Moveable left to right.**
- On multi joint arm for easy accessibility with holder for high and low vacuum suction, 3 way syringe and chair control.
- **BeeCool Plus:** Light cure cordless having 2modes, **Turbo mode:** 1s, 3s, illumination 2300 -2500mW/cm, **Normal mode:** 5s, 10s, 15s, illumination 1000-1200 mw, Constant light intensity. Metal light guide, With built in intensity meter.
- **Comfort Operator's Stool:** Zero Back Ache, Having height adjustment, **Backrest:** height adjustment and has swing movement, forward and backward movement, and also fixed position, supported by pneumatic piston. It moves with body giving support and giving certain pressure all the time on the back muscle which prevents professional hazard like hunch back/ back pain which is common in dentist, this has great role to play in long time clinical practice.

Right Hand Working



Left Hand Working

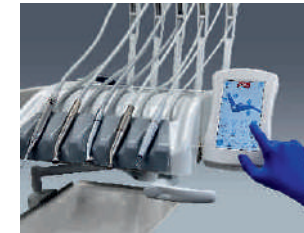

Marketed By:



MEDIKABAZAAR

Boston Ivy Healthcare Solutions Private Limited

3rd Floor, Second Avenue, Atul Projects,
Subhash Nagar, Andheri East, Mumbai,
Maharashtra 400093.

+91 794805 8625

support@medikabazaar.com


www.medikabazaar.com

CHAMUNDI TOUCH
DENTAL CHAIR MOUNT UNIT



Chamundi Touch Dental Chair Mount Unit

- **Programmable Dental Chair**, control with 7inch touch screen with HDMI all in one solution in chair, Powered by 64 bit quad core processor with 8GB RAM having 4 user control, Each user has 4 programmable chair position, Zero and gargle position.
- Control for water for spittoon and tumbler.
- **The Menu Screens:** Micromotor control, Scaler control, Airotor, Flush and Setting screen.
- **Memory Foam:** Seat, backrest and head rest are fitted with Most latest memory foam covered with Ultra Leather from USA.
- **Light:** Floor mount LED light ("FARO" MAIA) Operating light having On/Off & Intensity control by (un-touch) sensor.
- **Head Rest:** Having lateral left to right movement, rising and lowering, rotation on central axis with pneumatic rigid lock.
- **Hand Rest:** Euro moulded rigid, Dropping down for easy access of patient.
- **Spittoon Component:** Spittoon move automatically over the body, synchronised with gargle, and goes back to home position when the chair goes to working position. Also can be moved manually.
- Oval shaped bowl to cover more area for easy spitting, detachable spittoon, autoclavable, Special shape water component, Most graceful and clean.
- **Modular Trolley System:** Heavy duty trolley arm 90 degree shape with hinge with manual lock, well balanced.
- Over hanging delivery system, perfect locking by pull, stable, unlocking by pulling and leave the tubing.
- All Valves for air and water are solenoid controlled, Electrically operated 24 VDC valve for life long service.
- **NSK Brushless Optic Micromotor** - NLX Nano, fiber optic with internal spray having chipair function, supplied with NSK 1:5 speed increase contrangle handpiece having optic light and internal spray, push button type.
- **Airotor Points - 2Nos :**
 - a). Super torque Optic **NSK- Dyna LED** handpiece - 1 No
 - b). Mini head handpiece **NSK-Pico** smallest head in the world - 1 No
- 6 way syringe warm water and air - 1No. (LUZZANI Italy)
- 3 way syringe - 1No
- NSK Optic Scaler with 3 scaling tips - 1 Set
- **Multi-Function Foot Control**
- **Monitor Mounting Arm**
- **Wetline Suction:** continuous rating 300 L/Minute - 1No
- **Assistant Console:** on multi joint arm for easy accessibility with holder for high and low vacuum suction, 3 way syringe and chair control.
- **BeeCool Plus:** Light cure cordless having 2modes, **Turbo mode:** 1s, 3s, illumination 2300 -2500mW/cm, **Normal mode:** 5s, 10s, 15s, illumination 1000-1200 mw, Constant light intensity. Metal light guide, With built in intensity meter.
- **Handpiece Flushing System:** Used instruments can be flushed after surgery, two bottle with Separate switch for cleaning with disinfectant and clean water.
- **Suction Tube Cleaning:** We have designed an innovative system for cleaning suction tubes and filter to remove bio film.
- **Comfort Operator's Stool:** Zero Back Ache, Having height adjustment, **Backrest:** height adjustment and has swing movement, forward and backward movement, and also fixed position, supported by pneumatic piston. It moves with body giving support and giving certain pressure all the time on the back muscle which prevents professional hazard like hunch back/ back pain which is common in dentist, this has great role to play in long time clinical practice.

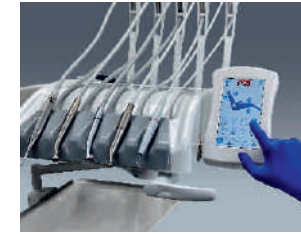

Marketed By:



MEDIKABAZAAR

Boston Ivy Healthcare Solutions Private Limited

3rd Floor, Second Avenue, Atul Projects,
Subhash Nagar, Andheri East, Mumbai,
Maharashtra 400093.

+91 794805 8625

support@medikabazaar.com


www.medikabazaar.com

CHAMUNDI PLUS
DENTAL CHAIR MOUNT UNIT



Chamundi Plus Dental Chair Mount Unit

- **Programmable Dental Chair** with 6 programme, 4 programme erasable, gargle, zero.
- **5" TFT based display** having 480/800 resolution, control panel screen with option for 4 programmable chair position, Zero and gargle position, Control for water for spittoon and tumbler.
- **The Menu Screens:** Micromotor control, Scaler control, Airotor, Flush and Setting screen.
- **Upholstery:** Seat, backrest and head rest are fitted with foam covered with rexine.
- **Light:** Floor mount LED light ("FARO" MAIA) Operating light having On/Off & Intensity control by (un-touch) sensor.
- **Head Rest:** Having lateral left to right movement, rising and lowering, rotation on central axis with pneumatic rigid lock.
- **Hand Rest:** Euro moulded rigid, Dropping down for easy access of patient.
- **Spittoon Component:** Spittoon move automatically over the body, synchronised with gargle, and goes back to home position when the chair goes to working position. Also can be moved manually.
- Oval shaped bowl to cover more area for easy spitting, detachable spittoon, autoclavable, Special shape water component, Most graceful and clean.
- **Modular Trolley System:** Heavy duty trolley arm U shape with hinge with manual lock, well balanced.
- Over hanging delivery system, perfect locking by pull, stable, unlocking by pulling and leave the tubing.
- All Valves for air and water are solenoid controlled, Electrically operated 24 VDC valve for life long service.
- **NSK Brushless Optic Micromotor** - NLX Nano, fiber optic with internal spray having chipair function, supplied with NSK 1:5 speed increase contrangle handpiece having optic light and internal spray, push button type.
- **Airotor Points - 2Nos :**
 - a). Super torque Optic **NSK- Dyna LED** handpiece - 1 No
 - b). Mini head handpiece **NSK-Pico** smallest head in the world - 1 No
- 6 way syringe warm water and air - 1No. (LUZZANI Italy)
- 3 way syringe - 1No
- NSK Optic Scaler with 3 scaling tips - 1 Set
- **Multi-Function Foot Control**
- **Monitor Mounting Arm**
- **Wetline Suction:** continuous rating 300 L/Minute - 1No
- **Assistant Console:** on multi joint arm for easy accessibility with holder for high and low vacuum suction, 3 way syringe and chair control.
- **BeeCool Plus:** Light cure cordless having 2modes, **Turbo mode:** 1s, 3s, illumination 2300 -2500mW/cm, **Normal mode:** 5s, 10s, 15s, illumination 1000-1200 mw, Constant light intensity. Metal light guide, With built in intensity meter.
- **Handpiece Flushing System:** Used instruments can be flushed after surgery, two bottle with Separate switch for cleaning with disinfectant and clean water.
- **Comfort Operator's Stool:** Zero Back Ache, Having height adjustment, **Backrest:** height adjustment and has swing movement, forward and backward movement, and also fixed position, supported by pneumatic piston. It moves with body giving support and giving certain pressure all the time on the back muscle which prevents professional hazard like hunch back/ back pain which is common in dentist, this has great role to play in long time clinical practice.

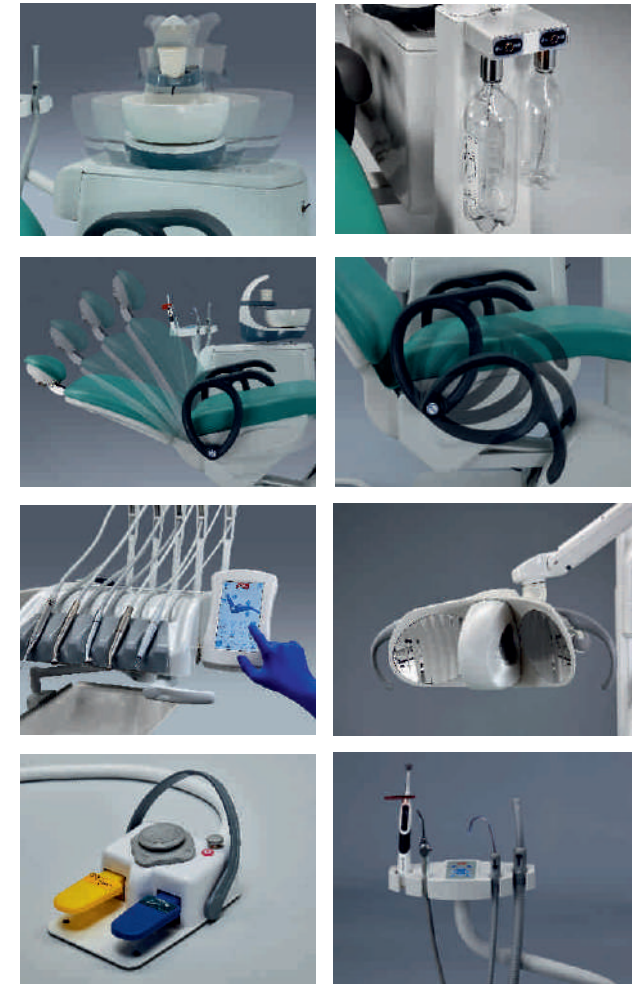

Marketed By:



MEDIKABAZAAR

Boston Ivy Healthcare Solutions Private Limited

3rd Floor, Second Avenue, Atul Projects,
Subhash Nagar, Andheri East, Mumbai,
Maharashtra 400093.

+91 794805 8625

support@medikabazaar.com

www.medikabazaar.com

CHAMUNDI BASIC
DENTAL CHAIR MOUNT UNIT



Chamundi Dental Chair Mount Unit

Body Contoured Electrically operated dental chair with state-of-art microprocessor based chair control system.

- Soft key switches for chair operation.
- 2 customizable programmes, zero position and gargle position.
- Synchronised movement of back rest and seat rest for patient comfort.
- Articulated Headrest with twin joint, Seat, backrest and head rest are fitted with foam covered with rexine.
- Swivelling arm rest for easy patient access.

Maia LED (Faro, Italy) Operating light having multi surface reflector having 35000lux.

- On / Off and intensity control by (un-touch) sensor with 3 directional movements.
- Light Fitted on well balanced base mount pillar & sturdy arm for operator's convenience.
- Monitor mounting arm

Easily cleanable chair side blue oval glass / round ceramic spittoon with auto water connection for spittoon and tumbler, manually moveable towards the patients chest position for convenience.

High and Low vacuum motorised non-stop CRS suction.

- Suction power of 300 Ltr / min.
- Auto drain & Auto flush system.

Assistant side Chair controls.

- High & Low vacuum suction cannulas with individual filters collect debris with removable tubing to keep hygiene.
- Beecool Plus LED Light cure unit (Cord less)
- 3-way syringe

Modular delivery system (Hanging / Over the Patient)

- Airtor points with tubing & nut fitting - 2 Nos
- Supreme Micro motor 35,000rpm with digital display of speed
- Supreme optic scaler with five scaling tips
- 3-way syringe
- SS Instrument Tray
- X-Ray viewer.

Multi-Function Foot Control

- Joystick function for chair movement, Airtor and Micromotor and scaler Operations

"COMFORT" Dental operator's stool moving on five castors with chromium plated base.

- Pneumatic piston for up down movement, backrest & seatrest.
- Ergonomically designed hand rest for operator's convenience.



capture the tiny details of the operative field with minimal effort, constantly keeping the colour temperature at 5,000 K.

FARO MAIA LED: Even the sunlight affects colour perception: at changing of the seasons or hours of the day we see different colours. Therefore it is important to have a light source that, at any time, is able to give accurate information to the doctor. All FARO lamps offer intensity adjustment so that the visual function can



Marketed By:

 **MEDIKABAZAAR**

Boston Ivy Healthcare Solutions Private Limited

3rd Floor, Second Avenue, Atul Projects,
Subhash Nagar, Andheri East, Mumbai,
Maharashtra 400093.

+91 794805 8625

support@medikabazaar.com

www.medikabazaar.com

CHAMUNDI BASIC
DENTAL CHAIR MOUNT UNIT



Chamundi Dental Chair Mount Unit

Body Contoured Electrically operated dental chair with state-of-art microprocessor based chair control system.

- Soft key switches for chair operation.
- 4 customizable programmes, zero position and gargle position.
- Synchronised movement of back rest and seat rest for patient comfort.
- Articulated Headrest with twin joint, Seat, backrest and head rest are fitted with foam covered rexine.
- Swivelling arm rest for easy patient access.

LED-C-20 Operating light having multi surface reflector with 50000lux.

- On / Off and intensity control by (un-touch) sensor with 3 directional movements.
- Light Fitted on well balanced base mount pillar & sturdy arm for operator's convenience.
- Monitor mounting arm

Easily cleanable chair side blue oval glass / round ceramic spittoon with auto water connection for spittoon and tumbler, manually moveable towards the patients chest position for convenience

High and Low vacuum motorised continues rating CRS suction.

- Suction power of 300 Ltr / min.
- Auto drain & Auto flush system.

Assistant side Chair controls.

- High & Low vacuum suction cannulas with individual filters collect debris with removable tubing to keep hygiene.
- **BeeCool Plus:** Light cure cordless having 2 modes, **Turbo mode:** 1s, 3s, illumination 2300 -2500mW/cm, **Normal mode:** 5s, 10s, 15s, illumination 1000-1200 mw, Constant light intensity. Metal light guide, With built in intensity metre.
- 3-way syringe

Modular delivery system over the patient

- Airtor points with tubing & nut fitting - 2 Nos
- a). Super torque Optic Dyna LED handpiece Supreme - 1 No.
- b). Miniature Handpiece MDK, Korea - 1 No
- **NSK Brushless Optic Micromotor** - NLX Nano, fiber optic with internal spray supplied with NSK 1:5 speed increase contrangle handpiece having optic light and internal spray, push button type.
- Supreme optic scaler with five scaling tips
- 3-way syringe
- SS Instrument Tray
- X-Ray viewer.

Multi-Function Foot Control

- Joystick function for chair movement, Airtor and Micromotor and scaler Operations
- **Comfort Operator's Stool:** Zero Back Ache, Having height adjustment, **Backrest:** height adjustment and has swing movement, forward and backward movement, and also fixed position, supported by pneumatic piston. It moves with body giving support and giving certain pressure all the time on the back muscle which prevents professional hazard like hunch back/ back pain which is common in dentist, this has great role to play in long time clinical practice.



Marketed By:



MEDIKABAZAAR

Boston Ivy Healthcare Solutions Private Limited

3rd Floor, Second Avenue, Atul Projects,
Subhash Nagar, Andheri East, Mumbai,
Maharashtra 400093.

+91 794805 8625

support@medikabazaar.com

www.medikabazaar.com

MEENAKSHI
DENTAL CHAIR MOUNT UNIT



Meenakshi Dental Chair Mount Unit

Body Contoured Electrically operated dental chair with state-of-art microprocessor based chair control system.

- Soft key switches for chair operation.
- 2 customizable programmes, zero and gargle programme.
- Articulated Headrest with twin joint, Seat, backrest and head rest are fitted with foam covered with rexine.
- Swivelling arm rest for easy patient access.

EDI LED (Faro, Italy) Operating light having multi surface reflector with capacity of 30000lux.

- On / Off and intensity control by (un-touch) sensor with 3 directional movements.
- Light Fitted on well balanced arm for operator's convenience.

Easily cleanable chair side round glass spittoon with auto water connection for spittoon and tumbler, Manually moveable towards the patients chest position for convenience.

High and Low vacuum motorised non-stop CRS suction.

- Suction power of 300 Ltr / min.
- Auto drain & Auto flush system.

Assistant side Chair controls.

- High & Low vacuum suction tubing with individual filters.
- Beecool Plus LED Light cure unit (Cord less)
- 3-way syringe

Modular delivery Hanging system.

- Airtor points with tubing & nut fitting - 2 Nos
- Supreme Micro motor 35,000rpm with digital display of speed
- Supreme optic scaler with five scaling tips
- 3-way syringe
- SS Instrument Tray
- X-Ray viewer.

Multi-Function Foot Control

- Joystick function for chair movement.
- Pedal operated airtor and micromotor control.

" EURO " Dental operator's stool moving on five castors with chromium plated base.

- Pneumatic piston for up down movement.
- Seat fitted with foam covered rexine.
- Backrest PU molded with height adjustment.



Marketed By:

 **MEDIKABAZAAR**

Boston Ivy Healthcare Solutions Private Limited

3rd Floor, Second Avenue, Atul Projects,
Subhash Nagar, Andheri East, Mumbai,
Maharashtra 400093.

 +91 794805 8625

 support@medikabazaar.com

 www.medikabazaar.com

MOOKAMBIKA
DENTAL CHAIR MOUNT UNIT



Mookambika Dental Chair Mount Unit

Body Contoured Electrically operated dental chair with state-of-art microprocessor based chair control system.

- Soft key switches for chair operation with Zero programme.
- Articulated Headrest with twin joint, Seat, backrest and head rest are fitted with foam covered with rexine.
- Swivelling arm rest for easy patient access.

LED C-20 Operating light having multi surface reflector with capacity of 50000lux.

- On / Off and intensity control by (un-touch) sensor with 3 directional movements.
- Reflector ensures the illumination uniformity of the pattern and good shadow less effect.
- Light Fitted on well balanced arm for operator's convenience.

Easily cleanable chair side round glass / porcelain spittoon with auto water connection.

- Auto drain & Auto flush system.
- High and Low vacuum motorised Dryco suction of 275Ltr / min.

Assistant side Chair console.

- High & Low vacuum suction tubing with individual filters.
- Beecool cord light cure unit
- 3-way syringe

Modular delivery Hanging system.

- Airotor points with tubing & nut fitting - 2 Nos
- Supreme Micro motor 35,000rpm with digital display of speed
- Supreme Non-optic scaler with five scaling tips
- 3-way syringe
- SS Instrument Tray
- X-Ray viewer.

Multi-Function Foot Control can use for Airotor and Micromotor Operating.

- Joystick function for chair movement.
- Pedal operated airotor & suction control.

" SURGI " Dental operator's stool moving on five castors with chromium plated base.

- Pneumatic piston for up down movement.
- Backrest & Seat fitted with foam covered rexine.



Marketed By:

 **MEDIKABAZAAR**

Boston Ivy Healthcare Solutions Private Limited

3rd Floor, Second Avenue, Atul Projects,
Subhash Nagar, Andheri East, Mumbai,
Maharashtra 400093.

+91 794805 8625

support@medikabazaar.com

www.medikabazaar.com

WAVE
DENTAL CHAIR MOUNT UNIT



Wave Dental Chair Mount Unit

Body Contoured Electrically operated dental chair with state-of-art microprocessor based chair control system.

- Soft key switches for chair operation.
- Articulated Headrest with twin joint, Seat, backrest and head rest are fitted with foam covered with rexine.

Light :

- Halogen bulb Operating light having glass reflector with two intensity
- On / Off and intensity control by switch.
- Light Fitted on Well balanced arm for operator's convenience.

Easily cleanable chair side round glass / porcelain spittoon.

Assistant side chair console

- Air ventury suction.
- Beecool cord light cure unit

Modular System

- Modular delivery Hanging system.
- Airtor point with tubing & nut fitting - 1No.
- Supreme Micro motor 35,000rpm with digital display of speed.
- Ultrasonic Non optic scaler with five scaling tip endo chuck - 1No.
- 3way syringe - 1No
- SS instrument Tray.
- X-Ray viewer.

Multi-Function Foot Control to operate Airtor & Micromotor.

Dental Operator Stool

- " SURGI " Dental operator's Stool moving on five castors with chromium plated base.
- Pneumatic piston for up down movement.



Marketed By:



MEDIKABAZAAR

Boston Ivy Healthcare Solutions Private Limited

3rd Floor, Second Avenue, Atul Projects,
Subhash Nagar, Andheri East, Mumbai,
Maharashtra 400093.

+91 794805 8625

support@medikabazaar.com

www.medikabazaar.com