



European Quality Parts

Work Mode Chart

	134° C Package	134° C Non Package	121° C Package	121° C Non Package	134B+
Sterilizing Temperature	134	134	121	121	134
Sterilizing Pressure	2.10	2.10	1.10	1.10	2.10
Sterilizing Time	6	4	20	16	18
Drying time min	15	15	15	15	15
Quick Drying time	1/5	1/5	1/5	1/5	1/5
Total time min	43 (29/33)	35 (21/25)	57 (43/47)	47 (33/37)	55 (41/45)

FENG 23L / SEA 23L-B-LED

- ④ B type fully automatic front-loading autoclave powered by triple vacuum cycle
- ④ Temperature Selection - 121° C / 134° C/ B+
- ④ Double door lock system for safety
- ④ 5 Programs for sterilization of all types of instruments – Solid/Hollow, Porous/Non-Porous etc.
- ④ Reduced full cycle time
- ④ Bow and Dick test qualified
- ④ User-friendly Control Panel
- ④ 2 Storage tanks for distilled and waste water
- ④ Data storage facility with USB functioning (2000 cycles)
- ④ Steam generator cleaning program to increase the efficiency of the autoclave
- ④ Chamber size 250 mm x 450 mm
- ④ High vacuum pump level -92 kPa
- ④ Optional micro printer

23 Ltr B-Type front loading Fully automatic with USB RECORDER




Aluminum Foil
Installation



High Vacuum
Pump



Stainless Steel
Connectors

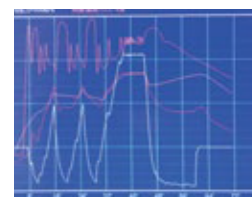

Reduced Full
Cycle Time